

Thank you for choosing our 4D Full Body Massager. Please read the instructions carefully before using the device.

This product is a household product and should not be used for commercial purposes. We will continually improve our products, so if there are any changes to the specifications, no further notice will be given.

PRODUCT OVERVIEW

The Probe 4D Full Body Massager combines advanced EMS micro-current, LED light therapy, and Radio Frequency Resurfacing to provide comprehensive skin and body care.

EMS Micro-Current

Widely used in the beauty industry, EMS micro-current massage helps accelerate blood circulation, effectively eliminate muscle fatigue and facial edema, and improve overall skin condition.

LED Color Light Therapy


Featuring three different wavelengths, the LED color lights can be used individually with the massage heads to enhance skin care. LED light therapy activates deep cells, accelerating skin metabolism, enhancing elasticity, and promoting firm, youthful skin.

Radio Frequency Resurfacing

This technology uses electromagnetic waves to stimulate subcutaneous collagen, causing it to shrink and regenerate over time. This results in firmer, wrinkle-free, and rejuvenated skin. The high-frequency vibration of water molecules in the cell layer generates heat (35°C-40°C / 95°F-104°F), causing collagen fibers to contract, which improves slack skin and reduces wrinkles. The thermal effect also promotes the regeneration and reorganization of collagen, repairing damaged and aging skin.

Experience the ultimate in skin and body care with the Probe 4D Full Body Massager, designed to deliver professional-grade results in the comfort of your home.

2

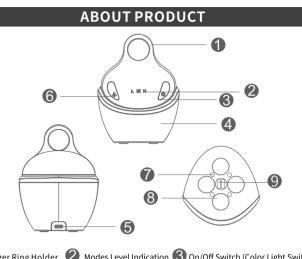


1 Finger Ring Holder 2 Modes Level Indication 3 On/Off Switch (Color Light Switch) **5** USB Socket 6 EMS/RF Switch

4 Charging Base

9 Charging Point

FUNCTION PARAMETER TABLE EMS (Micro-Current) 90Khz (L,M,H) RF (RADIO FREQUENCY) MODES 1MHz, Thermal Massage (L,M,H) LED Color Light Modes Blue Light, Yellow Light, Red Light, Flash Light, All Lights on High, Medium, Low Speeds Level

ELECTRICAL SCHEMATIC DIAGRAM

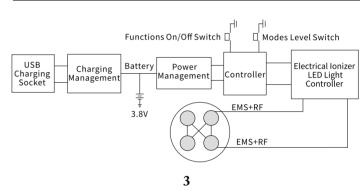

HOW TO USE


I:Short press the On/Off Switch, four red LED Lights turn on and the massage heads start working with low speed.









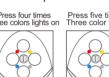



*Notes: When press the sixth time, LED lamp recycle from red light


II: Press the EMS/RF Switch, two switchable modes as shown

Massage Mode 4 massage heads simultaneous rotation Massage + EMS + RF Mode

L,M,H adjustable 4 massage heads simultaneous rotation Electrical Muscle Stimulation (EMS) High Radio Frequency (RF)

The massage speeds are adjustable with the Low speed, Medium speed, In Massage + EMS + RF Mode, the strength of EMS and RF increase along

Upon every press on EMS/RF Switch, Modes and strength level will change as below:

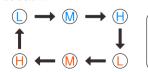

When letter "L, M, H" become blue color, Blue indicates Massage Mode; When letter "L" becomes red color, the Massage+EMS+RF Mode is active.

III: Body Parts Using Instruction

This product is recommended for use on the face, neck and shoulders. The methods can be referred as following picture.




Please apply facial gel to moisturize skin when selected the Massage + EMS

IV: Device will automatically shut down after 10 minutes of use.

Facial Care


- · Recommended usage time: 10 minutes per session, with 2-3 minutes in each area.
- Usage modes: Massage + LED mode, Massage + EMS + RF + LED mode.approximately 43°C / 109.4°F.
- **Note:** Please avoid using the device on the larynx. Since individual responses vary, adjust the massage probe to the appropriate facial spot during use.

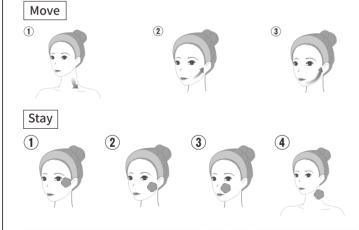

Neck & Shoulder Care

1

- · Recommended usage time: 10 minutes per session.
- Usage modes: Massage + LED mode, Massage + EMS + RF

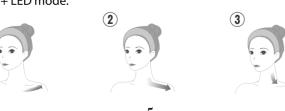

SPECIFICATION Product name Probe 4D

Probe 4D		
5V	Model No.	ET11
1050mAh	Rated Power	8W
227g	Charging time	5 hours
92×93×123mm	Working time	80 mins
	5V 1050mAh 227g	5V Model No. 1050mAh Rated Power 227g Charging time

Tips:

- 1. For best results, use the device with skincare products.
- 2. Ensure the massage heads stay in contact with your skin and adjust their positions to maximize effectiveness.

CAUTION

MAINTENANCE

- 1. Remove jewelry or ornaments before using the massager.
- 2. Use a dry towel dipped in water or a suitable cleansing liquid to clean the massage heads before and after use.
- 3. Charging Indicator (M): The light will be red while charging and turn green once fully charged.

SPECIAL REMINDER

- 1. Avoid pressing too hard on the skin, particularly in sensitive areas like the neck and chin.
- 2. Apply gel or skin care products to the skin before
- 3. 3. To prevent injury, do not insert your fingers between the massage heads while the device is in use.

6

OTHER INFORMATION

- Do not use or store this product in wet environments or at high temperatures.
- Tie up your hair before use to avoid it getting caught in the device.
- Keep the product out of reach of children.
- Do not disassemble or attempt to repair the product. Contact the appropriate after-service for assistance.

The following individuals should not use this product:

- People who are pregnant, breastfeeding, or menstruating.
- Those who have undergone plastic surgery or have fractures.
- Individuals with skin lesions, festering, or inflammation.
- Those with allergies to heat.
- People undergoing dental treatment.
- Individuals with sensitive skin.
- Those with hypertension, malignant tumors, blood disorders, internal diseases, tuberculosis, or an allergic constitution.

www.eternitech.shop info@eternitech.shop

7