
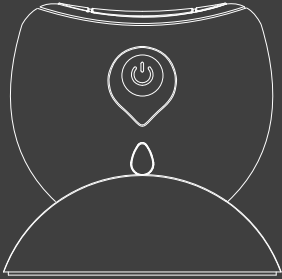
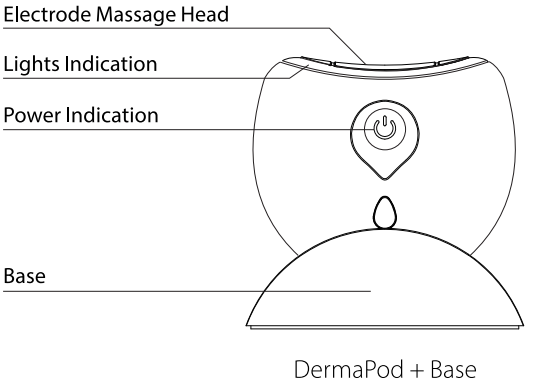
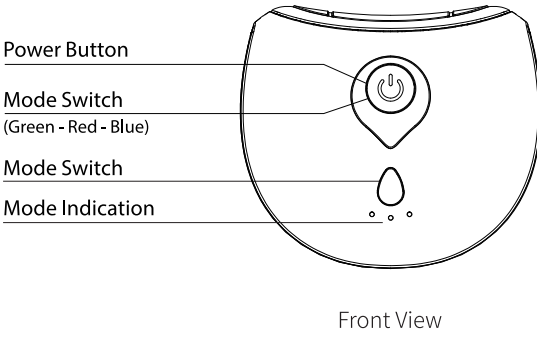
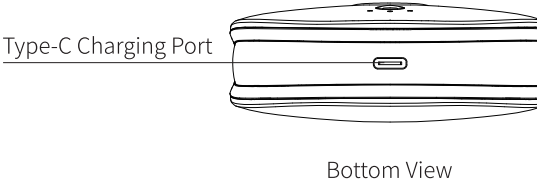
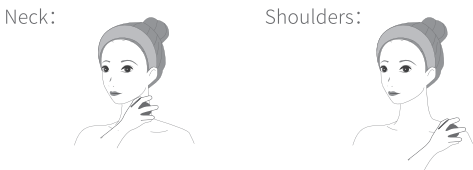
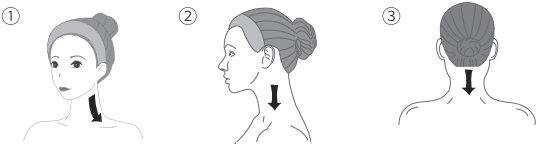
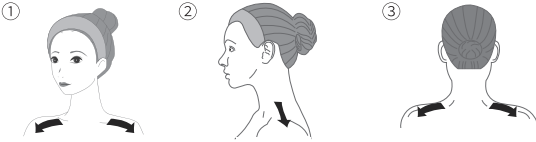


<div data-bbox="226 315 598 531"> ETERNITECH</div> <div data-bbox="296 670 528 902"></div> <div data-bbox="276 1020 550 1144"><p>DermaPod</p><p>Neck Therapy Device</p><p>User Manual</p></div>	<div data-bbox="742 226 986 264"><p>Product overview</p></div> <div data-bbox="742 294 1211 398"><p>Meet DermaPod, your new neck therapy device with a fashionable design that is easy to carry and operate.</p></div> <div data-bbox="742 441 1211 1023"><p>DermaPod uses EMS technology to reduce neck wrinkles, improve neck contours, and promote collagen synthesis. It can relieve neck muscle spasms, increase local blood circulation, and facilitate muscle retraining. This makes the skin more elastic, vibrant, and firm, improving sagging skin. The device combines EMS lifting and firming with constant temperature hot compress and high-frequency vibration massage, delivering a new beauty and skincare experience. DermaPod helps the skin better absorb nutrients, enhances muscle strength, increases blood circulation in the neck, and promotes local metabolism, thereby relieving neck muscle fatigue and eliminating soreness, stiffness, and other discomforts.</p></div>	<div data-bbox="1315 226 1562 264"><p>ABOUT PRODUCT</p></div> <div data-bbox="1315 277 1756 599"><p>Electrode Massage Head</p><p>Lights Indication</p><p>Power Indication</p><p>Base</p><p>DermaPod + Base</p></div> <div data-bbox="1315 675 1756 972"><p>Power Button</p><p>Mode Switch (Green - Red - Blue)</p><p>Mode Switch</p><p>Mode Indication</p><p>Front View</p></div> <div data-bbox="1315 1023 1756 1179"><p>Type-C Charging Port</p><p>Bottom View</p></div>	<div data-bbox="1859 226 2170 264"><p>Operation Instructions</p></div> <div data-bbox="1859 294 2287 617"><p>① Press power button once EMS + Massage mode, green light on.</p><p>② Press power button twice EMS + Massage + Hot Compress, red light on.</p><p>③ Press power button the 3rd time EMS + Massage + Cold Compress, blue light on.</p><p>④ Press power button the 4th time EMS + Massage + Cold Compress, blue light on.</p></div> <div data-bbox="1859 675 2337 743"><p>To turn off the device:</p><p>When the DermaPod is on, long press the power button.</p></div> <div data-bbox="1859 768 2270 861"><p>EMS Intensity Selection:</p><p>There are 3 EMS intensity levels. The red light will indicate the selected level.</p></div> <div data-bbox="1859 912 2337 1073"><p>Tips:</p><p>1. Cold Compress: The device takes about 20 seconds to reach approximately 19°C / 66.2°F.</p><p>2. Hot Compress: The device takes about 20 seconds to reach approximately 43°C / 109.4°F.</p></div>
<div data-bbox="184 1315 463 1353"><p>Gentle Care Routine</p></div> <div data-bbox="184 1378 630 1406"><p>It is recommended to use on neck and shoulders.</p></div> <div data-bbox="206 1449 595 1590"><p>Neck:</p><p>Shoulders:</p></div> <div data-bbox="184 1605 625 1668"><p>Please apply skin care products before having EMS + Massage + Cold / Hot compress.</p></div> <div data-bbox="184 1701 318 1728"><p>Neck Routine</p></div> <div data-bbox="194 1776 635 1897"></div> <div data-bbox="184 1942 368 1970"><p>Shoulders Routine</p></div> <div data-bbox="194 2016 635 2134"></div>	<div data-bbox="730 1315 1009 1353"><p>DermaPod Cautions</p></div> <div data-bbox="730 1378 1203 1713"><p>The following people should not use this product:</p><ul style="list-style-type: none"><li>• Pregnant or nursing individuals</li><li>• People with sensitive skin</li><li>• Individuals with skin lesions, festering wounds, or inflammation</li><li>• Those with heart pacemakers or other implanted medical devices</li><li>• Individuals with fever, heart disease, hypertension, malignant tumors, blood disorders, tuberculosis, or other serious conditions or mental health disorders. Please consult with a healthcare professional before use.</li></ul></div> <div data-bbox="730 1746 1203 1970"><p>Additional Precautions:</p><ul style="list-style-type: none"><li>• Keep the device away from children.</li><li>• Do not use diluents, benzene, solvents, or other chemicals to clean the product.</li><li>• Avoid dropping or striking the product. Stop using it if damaged.</li><li>• Do not disassemble or attempt to repair the product. Contact customer service for assistance.</li></ul></div> <div data-bbox="730 2003 1211 2255"><p>Maintenance:</p><ul style="list-style-type: none"><li>• Use a towel to clean any residue of skincare products left on the massager head.</li><li>• Recharge the product periodically if storing for an extended period.</li><li>• Store the product in a dry and ventilated place.</li><li>• Keep away from fire.</li><li>• Dispose of the product properly when it reaches the end of its lifespan.</li></ul></div>	<div data-bbox="1305 1736 1667 1819"><p>Website: <a href="http://www.eternitech.shop">www.eternitech.shop</a></p><p>Email: <a href="mailto:info@eternitech.shop">info@eternitech.shop</a></p></div>	<div data-bbox="1859 1315 2050 1353"><p>Specifications</p></div> <div data-bbox="1859 1383 2123 1680"><p>Product Name: DermaPod</p><p>Model: ET12</p><p>Rated Power: 8W</p><p>Rated Voltage: 5V</p><p>Battery Capacity: 1050mAh</p><p>Net Weight: 190g</p><p>Charging Time: ≤5H</p><p>Size: 90x88x54mm</p></div>